## NÁVOD NA POUŽITIE

Envy Therapy

Light Therapy  
Face Mask

RADICAL  
SKINCARE SYSTEM

ISKI

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MADE IN HONG KONG EXCLUSIVELY  
FOR ENVY THERAPY

DISTRIBUTOR:  
KLINIKA ENVY, S.R.O., SLOVAKIA

ENVYTHERAPY.COM

AT-HOME LED THERAPY IS A GREAT TOOL FOR IMPROVING THE FULL SPECTRUM OF SKIN CONDITIONS AND IMPERFECTIONS.

PHOTOTHERAPY IS USED BY DERMATOLOGISTS FOR ITS HIGH EFFICIENCY AND ITS RESULTS HAVE BEEN CONFIRMED IN NUMEROUS EXPERT STUDIES.

Envy Therapy

Light Therapy  
Face MaskRADICAL  
SKINCARE SYSTEM

|EN|

This instruction manual includes an explanation of the individual functions of the device as well as warnings when using it.

Before using the product, please read this instruction manual carefully.

Charge the device before using it for the first time.

This device is not a medical, but a cosmetic tool.

The effect of the device is influenced by the skin type and the type of cosmetic products used.

## TABLE OF CONTENTS

- PRODUCT INFORMATION
  - MAIN BENEFITS OF LED THERAPY
  - DESCRIPTION OF MODES
- USAGE
  - HOW TO USE LED MASK
  - OPERATING INSTRUCTIONS FOR THE REMOTE CONTROL
  - LED MASK EXPOSURE TIME
  - LED MASK FREQUENCY OF USE
- WARNINGS
- PACKAGE CONTENT
- PRODUCT SPECIFICATION

## 1. PRODUCT INFORMATION

## 1.1. MAIN BENEFITS OF LED THERAPY

At-home LED therapy is a great tool for improving the full spectrum of skin conditions and imperfections. Phototherapy is used by dermatologists for its high efficiency and its results have been confirmed in numerous expert studies.

It effectively eliminates signs of skin aging such as wrinkles and fine lines.

It reduces sebum production, acne, and is also suitable for sensitive skin.

Phototherapy helps with pigmentation, skin redness and uneven skin texture.


By stimulating collagen production, it makes the skin look firmer and younger.

It also improves skin texture and tone.

## 1.2. DESCRIPTION OF MODES

Envy Therapy LED Mask has eight light modes with different light colors. Each color has a different skin-improving clinically proven effect \*.

- Red light:** increases skin elasticity, stimulates cell metabolism, reduces wrinkles and fine lines.
- NIR – near-infrared light:** penetrates deeper into the tissue, promotes collagen and elastin production and tightens skin.
- Blue light:** decreases sebum production, eliminates bacteria causing acne, reduces inflammation and accelerates the healing process.
- Orange light:** corrects color irregularities of sun-damaged skin, deep pigmentation and melasma.
- Purple light:** soothes inflammatory processes, acts as an antioxidant, brightens the skin.
- Indigo light:** increases cellular energy, minimizes irritation of sensitive skin, speeds up post-treatment recovery.
- Green light:** reduces melanin production, evens skin tone, softens dark spots and other color imperfections.
- Yellow light:** positively affects the full range of rosacea symptoms and soothes skin redness and edema.

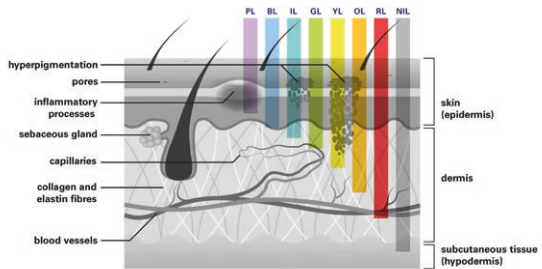PL (Purple Light) 400 nm  
BL (Blue Light) 465 nm  
IL (Turquoise Light) 495 nm  
GL (Green Light) 530 nm

YL (Yellow Light) 590 nm  
OL (Orange Light) 615 nm  
RL (Red Light) 640 nm  
NIL (NIR-Infrared Light) 830 nm

\* The Journal of Clinical and Aesthetic Dermatology

## 2. USAGE

## 2.1. HOW TO USE LED MASK


- Use the mask on cleansed and dry skin without make-up. Some ingredients in cosmetic products may interfere with the effectiveness of LED therapy.
- The package includes goggles for the eye area. Insert them from the inside of the mask. Put the eye goggle with "L" on the left eye hole. Make sure that the letter "L" is on the outer part of the left eye hole. Repeat the same step with the right ("R") eye goggle. However using the goggles is not required.
- Assemble fixing straps. Fixing straps do not distinguish between left and right. Make sure that you place the flannel side of the strap so that it faces inwards - towards your face. Attach the mask to your head with the straps.
- Connect the mask cable to the remote control.
- Press and hold the ON/OFF button on the remote control to turn the mask on and set the desired time interval of 5 or 10 minutes. The top MODE button sets the preferred LED light color.
- The mask will turn off automatically after the time interval you have set, or you can turn it off whenever you wish to by pressing the ON/OFF button.
- You can move around freely with the mask on.
- After use, wipe the inner silicone part of the mask with a dry or a damp cloth, dry it, and store it properly.
- Then continue with your daily skincare routine.

DO NOT USE THE MASK IF YOU BELONG TO ONE OF THE FOLLOWING GROUPS:

- Pregnant or breastfeeding women.
- Photosensitive people or people who are taking medication that can cause photosensitivity.
- People with genetic eye or skin diseases (such as lupus erythematosus, photosensitive eczema, albinism, etc.).

## 2.2. OPERATING INSTRUCTIONS FOR THE REMOTE CONTROL

- ON/OFF button / Time setting:** Press and hold the ON/OFF button to turn the device on/off. When the mask is on, the default mode is red light mode, the indicator RL is on, and the default time is 5 minutes. Press the button shortly to set the exposure interval to 5 or 10 minutes.
- MODE button:** is intended for switching between different lights. The default mode is red light mode, the indicator RL is on.
  - After pressing the MODE button once, the red light mode changes to the near-infrared light mode and the NIR indicator will switch on.
  - After pressing the MODE button for the second time, it changes the NIR light mode to the blue light mode and the BL indicator will switch on.
  - After the third press, the blue light mode changes to the orange light mode and the OL indicator switches on.
  - The fourth press changes the orange light mode to the purple light mode and switches on the PL indicator.
  - The fifth press changes the purple light mode to the indigo light mode and switches on the IL indicator.
  - After the sixth press, the indigo light mode changes to the green light mode and the GL indicator switches on.
  - After the seventh press, the green light mode changes to the yellow light mode and the YL indicator switches on.
- Battery power indicator:** When the remote control is charged and turned on, the corresponding mode indicator lights up steadily. Once the battery of the controller starts to lose power, the mode indicator will start flashing. The remote control, and the mask itself, will automatically turn off when the battery is empty.
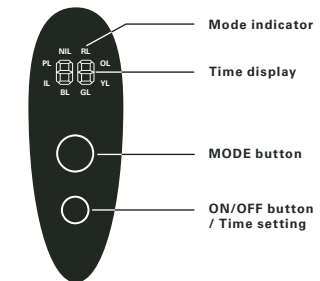
## 2.3. LED MASK EXPOSURE TIME

10 minutes.

## 2.4. LED MASK FREQUENCY OF USE



- We recommend using the mask 4 times a week in the first month.
- In the second month, use it 3 times a week.
- From the third month onwards, use the mask twice a week. You can continue to use the mask at this frequency indefinitely.

If you have multiple skin conditions, divide each treatment time into 5-minute intervals and choose your preferred light combination.

## 3. WARNINGS

- Please fully charge the remote control before first use and do not use it while charging.
- Do not clean, wash or use the remote control while charging. Avoid direct contact with a water source.
- If excessive redness or other adverse reactions occur, stop using the device immediately and consult a beautician or doctor.
- If you suspect that there is something wrong with the device, stop using it and contact the seller.
- Do not throw or thump the product.
- Do not store the mask in excessively high or low temperatures or in a damp environment. Store it in its original packaging in a cool, dry place. Do not expose it to direct sunlight.
- Do not place the device on blankets or electric blankets.
- Charge the device using only the original charging cable.
- Be careful when using and storing the mask. It is not suitable for use by children, keep it out of their reach.
- The device is not intended for consumption.
- Do not bring the product, remote control or the charging cable into contact with liquids or metal products.
- Unauthorized assembly, disassembly, repair or modification of the device may result in a malfunction that cannot be repaired and will void the warranty. Do not attempt to repair the device yourself and always contact the seller.

## 4. PACKAGE CONTENT
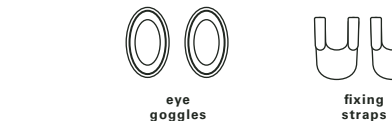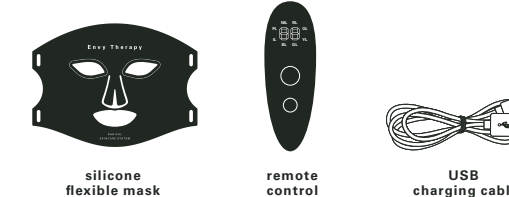
## 5. PRODUCT SPECIFICATION

PRODUCT MODEL: RB-008G  
NUMBER OF LEDs: 280 PCS  
CHARGING VOLTAGE: 5V  
NOMINAL VOLTAGE: 3,7V  
CHARGING CURRENT: 600MA  
NOMINAL POWER: 4W  
MADE OF FOOD GRADE SILICONE.

MADE IN HONG KONG EXCLUSIVELY FOR ENVY THERAPY

DISTRIBUTOR:  
KLINIKA ENVY, S.R.O., SLOVAKIA

ENVYTHERAPY.COM